


# Optilene<sup>®</sup> Mesh

Universal-lightweight mesh

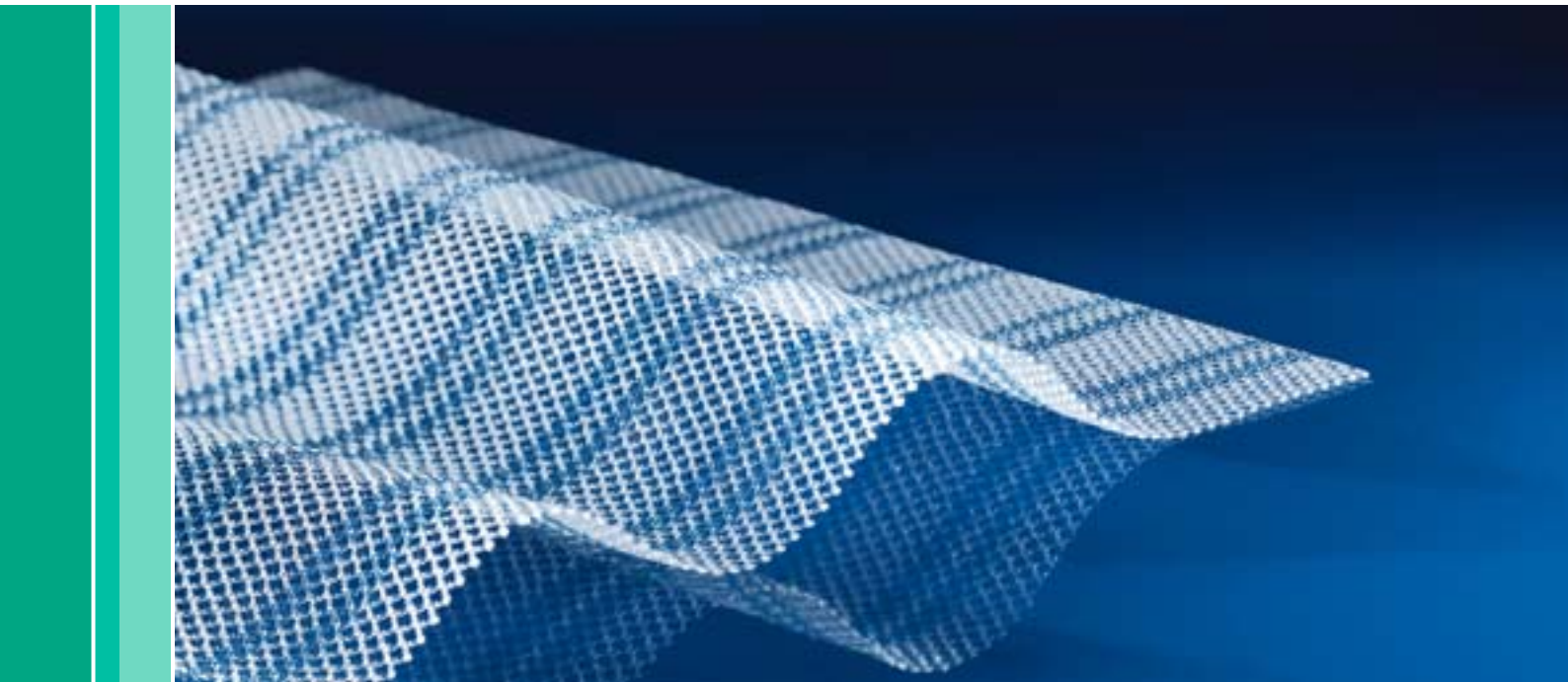

Hernia Repair

# Optilene® Mesh

Highest convenience for the patient and the surgeon

The new universal-light mesh combines a lightweight concept with an excellent handling. It is a soft, large pore and elastic mesh. At the same time it is a strong mesh with excellent spread characteristics.

When choosing Optilene® Mesh there is no need for changing the surgical technique. The innovative knit design of Optilene® Mesh provides optimal handling for open and laparoscopic application. The mesh can be used for inguinal and incisional hernia repair.



## universal-light

### ▶ **Universal-light mesh**

excellent handling lightweight mesh

### ▶ **Large pores**

for optimal healing and scar formation

### ▶ **Blue guidelines**

for a well directed placement

### ▶ **Filament-reinforced**

structure with ideal spread characteristics

### ▶ **Fully transparent**

to see the tissue below

### ▶ **Homogenous elasticity**

in longitudinal and latitudinal direction

### ▶ **Open and laparoscopic procedures**

one mesh suitable for all hernia repairs


# Large pores



## **Filament-reinforced**

The blue guidelines allow for a well directed placement of the mesh.

The lines are not only coloured but also designed to reinforce the mesh so that the mesh is able to selfexpand easily.

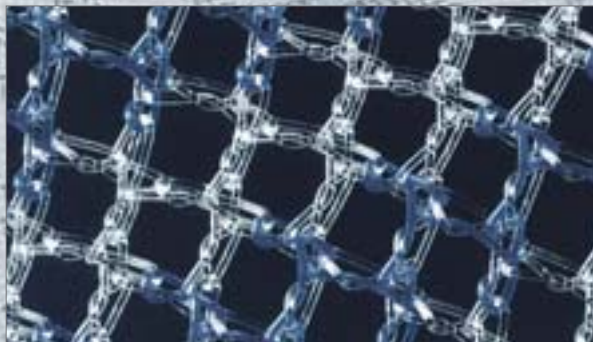

## **Pore size**

The large pore structure enables for optimal healing and formation of an elastic scar.

The innovative knit design offers full transparency and a homogenous elasticity.

## Product range

Article No.	Dimensions	Contents
1065020	5 x 10 cm	5 pieces
1065030	7.5 x 15 cm	5 pieces
1065040	10 x 15 cm	5 pieces
1065060	26 x 36 cm	5 pieces
1065080	15 x 15 cm	5 pieces
1065090	30 x 30 cm	5 pieces
1065140	4.5 x 10 cm	5 pieces
1065150	6 x 14 cm	5 pieces



# Optilene® Mesh


Highest convenience for the patient and the surgeon

## Spread characteristics

Overall flexural rigidity (Measurements according to DIN 53362)

The innovative structure of Optilene® Mesh combines tremendous weight reduction with excellent spread characteristics (weight reduction 30 - 40 % compared to traditional meshes).

Optilene® Mesh offers excellent handling be it in laparoscopic or open hernia repair.

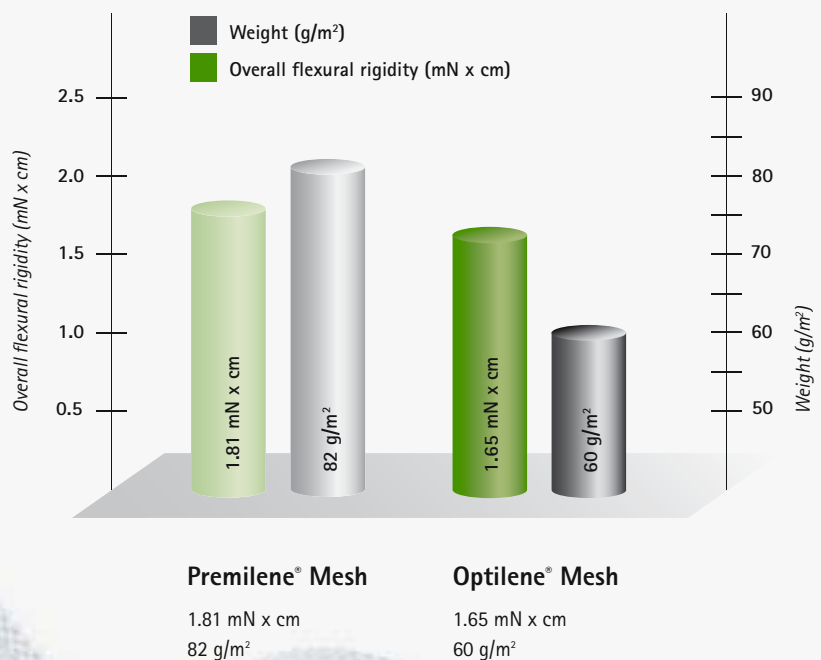

## Product features

<b>Material</b>	Monofilament Polypropylene
<b>Construction</b>	Knitted
<b>Weight</b>	60 g/m <sup>2</sup>
<b>Pore size</b>	1.5 mm
<b>Thickness</b>	0.53 mm
<b>Indications</b>	Inguinal Hernia Incisional Hernia Reconstruction of chest wall
<b>Contraindications</b>	Optilene® Mesh should not be implanted in contaminated and infected areas as well as in children during the growth phase.
<b>Sterilisation</b>	Ethylene oxide

Optilene® Mesh



B. Braun Surgical SA | Carretera de Terrassa, 121 | 08191 Rubí | Spain  
Phone +34 93 5 86 62 00 | Fax +34 93 6 99 73 03 | [www.bbraun.es](http://www.bbraun.es)

Aesculap AG | Am Aesculap-Platz | 78532 Tuttlingen | Germany  
Phone +49 (0) 74 61 95-0 | Fax +49 (0) 74 61 95-26 00 | [www.aesculap.com](http://www.aesculap.com)

Aesculap – a B. Braun company

All rights reserved. Technical alterations are possible. This leaflet may be used in no way whatsoever other than for the offering, buying and selling of our products.  
No part may be copied or reproduced in any form. In the case of misuse we retain the right to recall our catalogues and pricelists and to take legal action.