



# CST STAPLER

INTRALUMINAL STAPLER FOR SINGLE USE

# CST STAPLER

INTRALUMINAL STAPLER FOR SINGLE USE



# PRODUCT FEATURES AND BENEFITS



## AN AUDIBLE BELL WHEN FIRED


Clear and distinctive via an audible signal.


## SMART ASL™ AUTOMATIC SAFETY LOCK

Innovative design engages and disengages safety lock automatically, reducing the risk of accidental misfiring and refiring.

## PRS™ PRESSURE RELEASE SYSTEM

Reduce tissue tension due to increased housing capacity, helping to eliminate the over compression of tissue, the "piston effect".

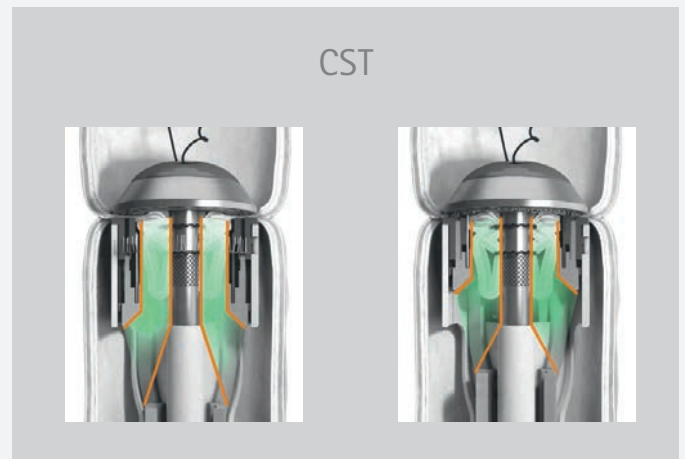

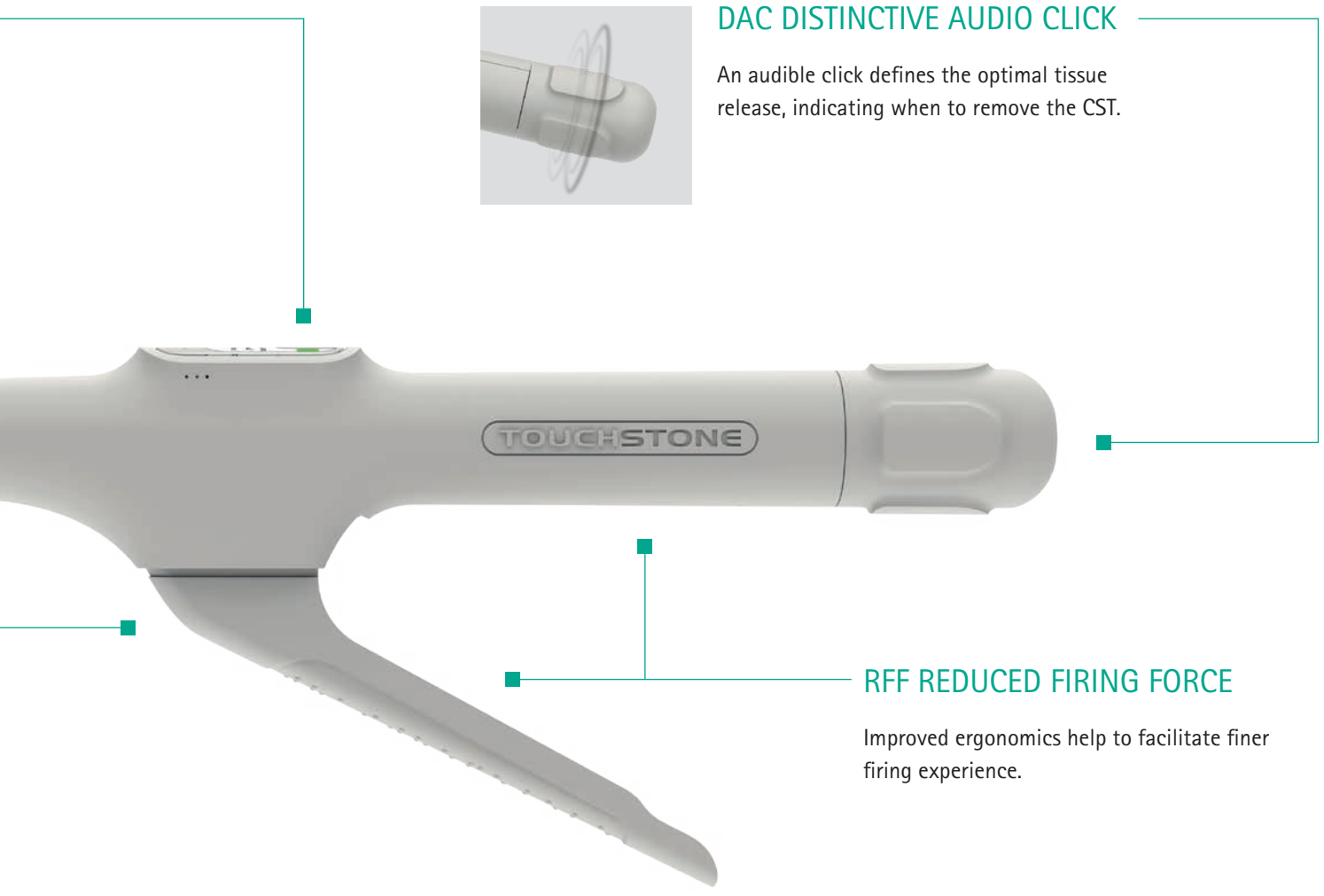


### DAC DISTINCTIVE AUDIO CLICK

An audible click defines the optimal tissue release, indicating when to remove the CST.

### RFF REDUCED FIRING FORCE

Improved ergonomics help to facilitate finer firing experience.

### OTHER

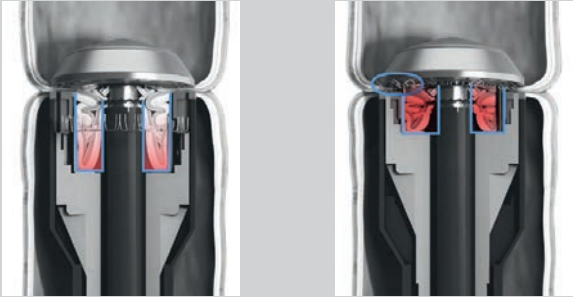

## ORDERING INFORMATION



PRODUCT CODE	CST21	CST25	CST29	CST33
DESCRIPTION	Barrier Free, high volume housing, air sealed circular stapler			
COLOR CODE	Orange	White	Blue	Green
STAPLE QUANTITY	16	20	24	28
ANVIL DIAMETER	21 mm	25 mm	29 mm	33 mm
BLADE DIAMETER	13.9 mm	17.0 mm	20.5 mm	24.8 mm
HOUSING VOLUME BEFORE/AFTER FIRING	5.1/4.6 ml	9/8.2 ml	12/10.9 ml	14.6/13.2 ml

B. Braun Surgical, S.A.U. | Carretera de Terrassa, 121 | 08191 Rubí | Spain  
Tel. +34 935 86 62 00 | +34 935 88 10 96 | [www.bbraun.es](http://www.bbraun.es)

Manufacturer acc. to MDD 93/42/EEC:  
Touchstone International Medical Science Co., Ltd. | 278 Dongping Street | Suzhou Industrial Park | 215123 Suzhou | China